




## Tracheal tube for long-term ventilation

- Unique *SealGuard™* high volume / low pressure cuff helps to ensure an efficient low pressure cuff seal with significant lower risk of aspiration.
- Specially moulded hooded tip to assist intubation and provides high patient safety and comfort.
- Murphy eye incorporated as an additional safety feature.
- Suitable for both oral and nasal intubation.

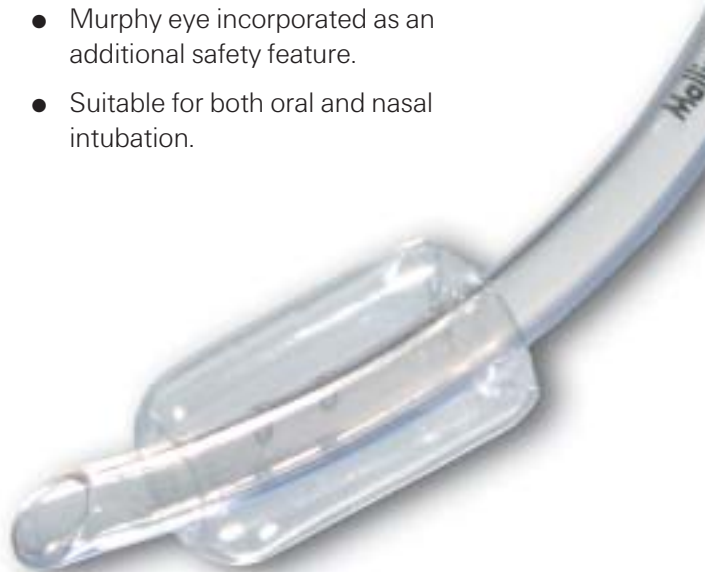

Due to the superior sealing properties of the unique cuff material a large amount of secretions could potentially build up above the cuff and possibly flow in to the pharynx. To minimize the potential risk of infection from the build up of secretions in the subglottic space, the *SealGuard™ Evac* tube is recommended. The *SealGuard™ Evac* tube incorporates an additional lumen for improved access and drainage of the subglottic space.

REF	I.D. mm	O.D. mm	Cuff Resting Diameter mm	Length mm
109650	5.0	6.9	18	274
109655	5.5	7.5	18	304
109660	6.0	8.2	21	317
109665	6.5	8.8	21	327
109670	7.0	9.6	24	339
109675	7.5	10.2	26	350
109680	8.0	10.9	28	360
109685	8.5	11.5	30	370
109690	9.0	12.1	32	375
109610	10.0	13.5	34	375

Trademarks of Tyco Healthcare Group LP or an affiliate © 2004 Tyco Healthcare Group LP. All rights reserved.

C-AN-SGuard-H1/GB

**tyco**  
Healthcare

MALLINCKRODT

